


## SHARP INNOVATION INJECTION NEEDLE



A simple, elegant solution, the 5 mm long lancet Injection Needle is offered in 22 and 25 gauge. Designed for single-handed operation, the innovative design limits accidental needle deployment during insertion and removal. The patented ratchet mechanism provides audible confirmation of deployment and minimizes needle bounce. A quick release button returns the needle fully into the sheath. The smooth metal hub and firm catheter aid insertion and ensures needle performance through various scope positions.

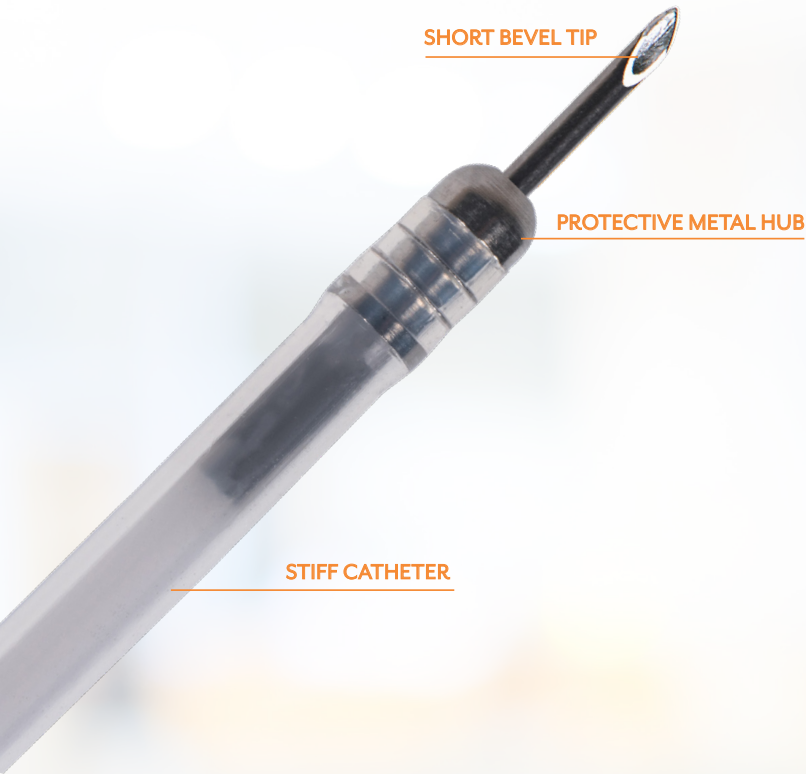

NEVER  
COMPROMISE  
ON QUALITY

DRAMATICALLY  
IMPROVE YOUR  
BOTTOM LINE

THROW AWAY  
CONTRACTS  
FOREVER

RELIABLE  
SUPPLY  
PARTNERSHIP

## KEY BENEFITS

### ONE-HANDED OPERATION

Ergonomically designed one-handed operation enables confident needle deployment that enables rapid needle retraction.

### PATENTED RATCHET DESIGN

Provides audible confirmation and secure locking of the deployed needle during use. As the needle is advanced, the ratchet gives tactile feedback and helps to reduce needle bounce.

### SHORT BEVEL TIP

The needle features a shorter 3-bevel tip, designed to improve accuracy when injecting submucosally. The smooth metal hub and firm catheter may support use for tamponade.

## SPECIFICATIONS

### INJECTION NEEDLES

Order Number	Henry Schein Item Number	Description	Needles Gauge/Length (mm)	Sheath Diameter (mm)	Working Length (cm)	Minimum Channel Size (mm)	Package Units
IN33131	129-0169	Injection Needle Upper	22/5	2.3	180	2.8	10/box
IN33221	129-0171	Injection Needle Upper	25/5	2.3	180	2.8	10/box
IN33151	129-0170	Injection Needle Lower	22/5	2.3	230	2.8	10/box
IN33241	129-0172	Injection Needle Lower	25/5	2.3	230	2.8	10/box

## CAN'T FIND IT?

Additional items may be available. Contact us if you can't find what you need.



**EXCEPTIONAL**  
ENDOSCOPIC  
SOLUTIONS.

2855 Boardwalk Drive  
Ann Arbor, MI 48103  
877.552.4027

**MTENDOSCOPY.COM**

CustomerCare@MTEndoscopy.com  
MPW33000 Rev. 3